

NEW

# Neuro Physical Training™

From Rehabilitation to Health Sport



# Neuro Physical Training™


## EXCLUSIVE HUBER® 360 TECHNOLOGY

HUBER® 360 features the new Multiaxis Motorized Platform. It comes equipped with force sensors built into the platform and handles. Trained by a health professional, your patients will work safely on both physical and cognitive skills: that's Neuro Physical Training! Once mobility and flexibility have been recovered on HUBER® 360, patients can enroll into a Health Sport program in your center to maintain and improve their health.

### MULTIAXIS MOTORIZED PLATFORM™ WITH BUILT-IN FORCE SENSORS

- ▶ INTEGRATED FUNCTIONAL ASSESSMENT ADAPTED TO ALL PATIENT TYPES
- ▶ DYNAMIC POSTURE CORRECTOR™ FOR A PRECISE AND PROGRESSIVE WORKOUT
- ▶ MULTI-DIRECTIONAL STIMULATION AND TARGETED RECRUITMENT

### SPECIFIC TRAJECTORIES FOR EACH THERAPEUTIC INDICATION

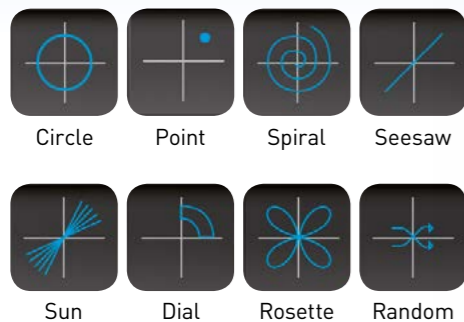

## 360° REHABILITATION

With HUBER® 360, the treatment is broken down into four movement fundamentals: Flexibility and Mobility, Dynamic Strengthening, Posture and Balance, and Resistance

### FLEXIBILITY AND MOBILITY

- ▶ The Dynamic Posture Corrector™ improves stretching precision.
- ▶ The Multiaxis Motorized Platform™ mobilizes joints in every plane of motion.



### DYNAMIC STRENGTHENING

- ▶ The Dynamic Posture Corrector™ insures proper positioning during exercise.
- ▶ Feedback on the screen allows control of strength and direction for better recruitment.



### POSTURE AND BALANCE

- ▶ The static and dynamic force platform improves objectification of posture.
- ▶ Balancing games are both fun and challenging for the patient.



### RESISTANCE

- ▶ Dynamic exercises adapted to the patient's physical condition increase their resistance.
- ▶ Short-duration high-intensity activities for faster results.



**NEW**

# HUBER® 360 TABLET APPLICATION

**A real patient follow-up tool adapted to your practice.  
Guaranteed to be an easy-to-use time saver.**

- › Full built-in assessment
- › Photo and video analysis tool
- › Connect to the machine wirelessly and without a USB drive
- › Send the assessment directly as a PDF to patients or prescribers
- › Create your own programs and free menus
- › HUBER® 360 Remote control
- › Monitor patients' exercise remotely



#### SIZE :

L=1800mm  
W= 1050mm (without guardrails)  
W=1335mm (with guardrails)  
H=2080mm

**WEIGHT :** 274 kg (with guardrails)

#### FLOOR SIZE :

1,90m<sup>2</sup> (without guardrails)  
2,4 m<sup>2</sup> (with guardrails)

#### POWER REQUIREMENTS :

230V / 50-60 Hz / 2 kVA

**MAXIMUM PLATFORM TILT :** 10°

#### MAXIMUM PLATFORM SPEED :

1 rev/s at 10° tilt

**MAXIMUM PLATFORM LOAD :** 150 kg

10.4" color screen



**ARTICULATED  
STOOL**



**PADS**



**REPOSITIONING  
BLOCK**



**HEART RATE MONITOR**

**LPG SYSTEMS S.A.**  
**INTERNATIONAL // MARKETING**  
2753, Route des Dolines - BP 243  
06905 Sophia-Antipolis Cedex France  
Tel.: +33 (0)4 92 38 39 00  
Fax: +33 (0)4 92 96 09 65  
[www.huber360.com](http://www.huber360.com)

